

Intelligent Evaporation



Hei-VAP Rotary Evaporators

Leading Safety Standards

Superior Ease of Use

Reduced Cost of Ownership






Do you demand the highest safety standards, perfectly engineered ease of use and reduced cost of ownership?

Intelligent Evaporation...

... stands for an overheat protection design that prevents water baths running dry, provides a detachable panel for superior ease of use from outside closed fume hoods and offers a proven vacuum seal and vapor tube system guaranteed for years of operation!

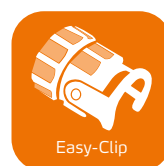

Leading **Safety Standards**

- Two independent safety circuits **avoid an overheat situation** of your heating bath: If bath temperature overshoots by 5 °C and if your water bath runs dry, heating will be switched off automatically
- **To avoid** an overtemperature situation and potential **thermal damage** to your sample, you can individually pre-program a maximum temperature on the Hei-VAP Precision, at which the bath powers off
- In case of a power cut, all motor lift models will remove the flask from the heating bath **preventing safety hazards** and potential thermal damage to your sample
- **To protect you from splashes and scalding**, the Hei-VAP heating bath features an integrated pour spout to remove bath fluid safely
- Heating bath is constructed with a **double wall** for your protection against burns and scalding
- The **guard hood covers the bath** even if the heating bath is adjusted horizontally to its most extreme limit
- A **metal guide** between the heating bath and the base unit **prevents bath instability**
- To **protect** your heating bath **against short-circuits and systematic corrosion**, the cable coupling complies with the protection class IP 67
- For your protection, the operating panel is located underneath a protective canopy so as to **keep your hands out of the dangerous zone** of rotating flasks, steam and splashing bath fluids
- All heating baths feature **non-slip handles** with safety grip to protect you from scalding bath liquids and allow for easy installation/removal of the bath
- A separate on/off switch for heating **prevents unintentional heat-up**. The button is illuminated for visual control
- In case of an overshoot of any parameter setting, all models will show a **warning on the operating panel**
- Glassware assemblies are optionally available with **transparent plastic safety coating** ensuring ultimate protection while allowing visual control at all times
- **IQ-OQ documentation** is available for Hei-VAP models on request

➤ Hei-VAP heating bath



YOUR ADVANTAGES

- The Hei-VAP heating bath reaches a temperature of 210 °C
- The universal heating bath can even accommodate 5-liter flasks
- The heating bath can be adjusted horizontally up to 200 mm
- To prevent short circuits and corrosion, the cable coupling complies with the protection class IP 67
- Non-slip handles for safe transportation eliminate the risk of scalding
- An integrated pour spout to remove bath fluid safely

Superior **Ease of Use**

- The detachable operating panel enables convenient handling from **outside closed fume hoods**
- The patented evaporating flask **clamp Easy-Clip cannot be lost** due to a connection to the vapor tube coupling ring
- The integrated screwed ring allows you to **remove sticking evaporating flasks** from the vapor tube **just in seconds** in a safe manner
- The universal heating bath can be adjusted horizontally up to 200 mm and offers efficient space to **accommodate even 5-liter flasks**
- The innovative ventilation and replenishment valves no longer have ground joints, which **enables maximum tightness, grease-free operation and prevents glass breakage**
- For your convenience, **the height of 160 mm and the angle of the evaporating flask can be adjusted easily from 20° to 80°**
- The optional automatic module Distimatic Bench-top allows for **continuous and unattended automatic operations**
- All components in direct or indirect contact with the medium are equipped with **FDA-compliant materials for an extended application range**






Reduced Cost of Ownership

- Achieve up to **40 % higher evaporation rates** due to unique tightness and industry-leading condensing surface with the **XL condenser**
- The high resistant and durable **PTFE vacuum seal achieves maximum tightness** and reduces spare part budgets
- The patented clamping sleeve ensures a **long life performance of the vapor tube**
- Increase your range of applications by using 5-liter flasks in the universal bath at **no additional cost**
- Maintenance-free and non-sparking motors **reduce repairs and down times significantly** to ensure years of continuous operation
- The flange on the condenser side of the drive is made of PPS (Polyphenylene sulfide), a robust industrial-designed resin with **superior chemical resistance** which provides more chemical resistance than stainless steel or welded aluminum. This reduces corrosion to a minimum which results in additional savings on maintenance and replacement components
- Spend your working hours more effectively on other challenges which require your attention and **evaporate large quantities continuously** with the automatic module Distimatic Bench-top
- Maximize your initial investment significantly: the automatic module Distimatic Bench-top is an **affordable and more flexible alternative** to large-scale rotary evaporators
- Reduce your maintenance costs. The sealed housing protects your rotary evaporator from aggressive fumes, liquids and vapors to prevent internal corrosion. This results in an **increased lifespan of 10 years** on average at a reduced maintenance and repair cost

➤ Vapor tube with patented clamping sleeve

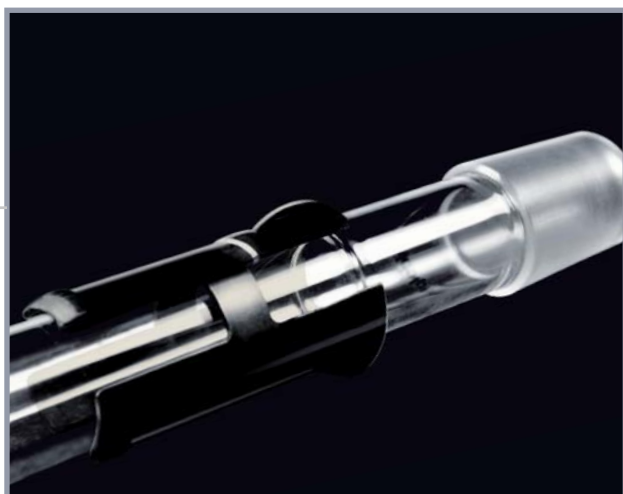

YOUR ADVANTAGES

- The patented clamping sleeve system allows easy removal of the vapor tube from the drive
- Years of performance life reduce your spare parts budget up to 75 %
- The Hei-VAP vapor tube does not get stuck in the drive, thus eliminating the risk of broken glass
- Immediate removal of the vapor tube from the drive facilitates your work

➤ Patented flask clamp Easy-Clip



YOUR ADVANTAGES

- Easy-Clip can't get lost due to a direct connection to the drive
- The integrated screw thread allows an easy removal of sticking flasks
- The Easy-Clip is made of highly resistant synthetic material which ensures a long life performance

➤ Operating panel



YOUR ADVANTAGES

- Changing all parameters from outside closed fume hoods guarantees the highest safety level
- Standard for all models is a panel with 1 m cable connection. Thus data security is guaranteed at any time
- Immediate access to the panel ensures ultimate control of the process while increasing your ease of use
- A separate on/off switch for heating prevents unintentional heat-up – the button is illuminated for visual control

Intelligent Evaporation



Leading Safety Standards

Superior Ease of Use

Reduced Cost of Ownership

The average operational **lifespan of 10 years** is backed by a **3 year warranty** and makes your purchase a worthwhile investment.

New, FDA-compliant, transparent and uv-resistant tubing, with **higher chemical resistance**

Achieve **up to 40 % higher evaporation rates** due to tightness and industry-leading **2,200 cm²** condensing surface **with the XL condenser**

The new highly resistant PTFE vacuum seals allow for long lasting constant application and **reduce spare parts budgets**. Also available as FFKM version for application with extremely aggressive solutions

Superior ease of use due to the patented Easy-Clip, which provides an **easy removal of sticking flasks** and non-sticking vapor tube with patented clamping sleeve that **eliminates broken glass** and reduces spare parts cost significantly

The innovative ventilation and replenishment valves no longer have ground joints, which enables **maximum tightness, grease-free operation and prevents glass breakage**

The detachable operating panel enables **safe and easy handling** from outside closed fume hoods

210 °C Safety heating bath with handles and pour spout will power off if temperature overshoots by 5 °C or runs dry

Bath power cable coupling complies with protection class IP 67 to **prevent corrosion and short circuits**

For unlimited evaporation:

In conjunction with the automatic module Distimatic Bench-top allows you unlimited and automated distilling in 24/7 operation. An affordable alternative to large-scale rotary evaporators

A separate on/off switch for heating **prevents unintentional heat-up** – the button is illuminated for visual control

A metal guide between the heating bath and the base unit **prevents bath instability**



➤ Hei-VAP Value / Hei-VAP Value Digital

Superior ease of use due to direct control

The Hei-VAP devices with hand lift are optimally suited for all standard applications. Controlling parameters is done so with the simple turn of a knob

- The convenience of large dial controls for adjustment of rotation speed and bath temperature
- **Hei-VAP Value Digital** models come with a digital display for set and actual heating bath temperature



Hei-VAP Value Digital HL/G3
P/N 560-01302-00

Model	Hei-VAP Value hand lift (HL)		Hei-VAP Value Digital hand lift (HL)	
	standard	coated	standard	coated
G1 Diagonal	560-01100-00	560-01110-00	560-01102-00	560-01112-00
G3 Vertical	560-01300-00	560-01310-00	560-01302-00	560-01312-00
G3 XL Vertical	560-01305-00	560-01315-00	560-01307-00	560-01317-00
G5 Dry Ice Condenser	560-01500-00	560-01510-00	560-01502-00	560-01512-00
G6 Vertical for Reflux	560-01600-00	560-01610-00	560-01602-00	560-01612-00

➤ Hei-VAP Advantage

Wide range of possibilities for reproducible results

With the Hei-VAP Advantage you have a direct view on all parameter settings and are able to expand your application range



Large digital LCD display which allows you to monitor all parameters even from a distance

Convenient motor lift adjusts height by the press of a button

- Available in hand lift or motor lift configurations
- Experience the advantages of our large digital 3.5" LCD display which indicates rotation speed, heating bath and vapor temperature
- Experience the advantages of the automatic process timer which turns off your evaporator at a pre-programmed time and removes evaporating flasks from bath on motor lift models
- Upgrade your Hei-VAP Advantage to a Hei-VAP Precision:

**Upgrade-KIT
Advantage Precision-HL**
P/N 569-30009-00

**Upgrade-KIT
Advantage Precision-ML**
P/N 569-40009-00

Hei-VAP Advantage ML/G3 XL
P/N 562-01305-00

Model	Hei-VAP Advantage hand lift (HL)		Hei-VAP Advantage motor lift (ML)	
	standard	coated	standard	coated
G1 Diagonal	561-01100-00	561-01110-00	562-01100-00	562-01110-00
G3 Vertical	561-01300-00	561-01310-00	562-01300-00	562-01310-00
G3 XL Vertical	561-01305-00	561-01315-00	562-01305-00	562-01315-00
G5 Dry Ice Condenser	561-01500-00	561-01510-00	562-01500-00	562-01510-00
G6 Vertical for Reflux	561-01600-00	561-01610-00	562-01600-00	562-01610-00

➤ Hei-VAP Precision

Maximum capabilities due to integrated vacuum control

The Hei-VAP Precision provides you with the maximum range of parameters, including but not limited to, reproducible applications and fully automated processes. Even solutions with different boiling points can be identified and distilled

- Available in hand lift or motor lift configurations
- Experience the advantages of our large 4.3" color graphic display which features all parameters, integrated vacuum controller and includes programs for automated distillations
- Make your life easier by utilizing the integrated software support which leads you through the capabilities of the evaporator
- For data management Hei-VAP Precision models feature a USB interface
- Save time by configuring the unit - all connected devices will be recognized by the evaporator, which eliminates the need to configure the unit manually
- Now you can manage even the most demanding applications easily and prevent foaming, excessive bubbling or bumping
- Save time with your individualized pre-programmed parameter settings for your 30 most common applications used on a daily basis by saving them in a "Favorites" section
- Each ramp offers 20 separate steps for demanding applications
- Apply an individual name to those ramps to avoid a mix-up
- Experience the advantages of the automatic process timer which turns off your evaporator at a pre-programmed time by raising the flask on motor lift models and releasing vacuum
- Control your process individually by pre-programming vacuum gradients
- Hei-VAP Precision models feature a built-in vacuum controller as a standard

For increased safety all bath temperature adjustments can be limited

Even from a distance you will be able to monitor the numbers on the display easily


Large color graphic display for all parameters and automated programs



Hei-VAP Precision ML/G3
P/N 564-01300-00

➤ Reduce your process times significantly via automated distillation programs


For automated distillation choose between these three programm opportunities:

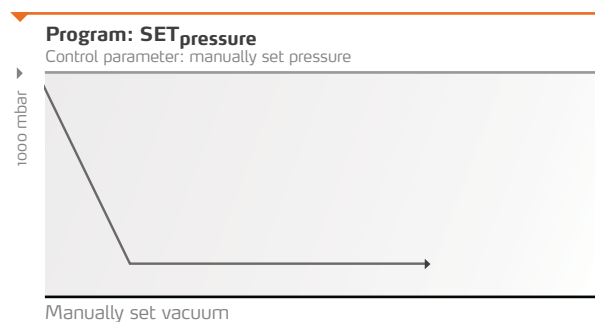

SET_{pressure}

Program for solvent separations


This program holds a pre-set vacuum constant. During the evacuation you have the full flexibility to choose another level of vacuum by the press of a button. Optionally, the process can be predefined via programmable ramps

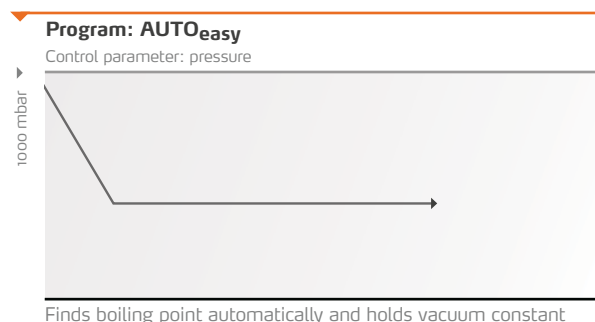

AUTO_{easy}

Supports the automatic pressure control for solvent separations

The pressure level for the first boiling point is found automatically. This program can be used with all Heidolph glassware assemblies and does not require any accessories

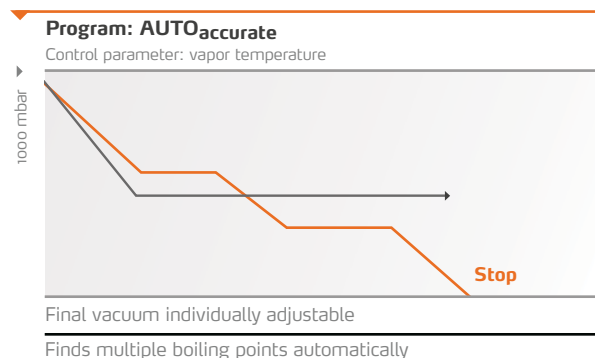

AUTO_{accurate}

Supports very precise automatic pressure control for solvent separations

This program works similar to AUTO_{easy} but in addition is able to identify multiple boiling points automatically by reducing vacuum accordingly. Also the AUTO_{accurate} program delivers the finest and most precise vacuum control. To take advantage of this program, the optional AUTO_{accurate}-sensor is required. This program cannot be used on Heidolph G1 and G5 glassware assemblies. Leading parameter of this program is the vapor temperature

Model	Hei-VAP Precision hand lift (HL)		Hei-VAP Precision motor lift (ML)	
	standard	coated	standard	coated
G1 Diagonal	563-01100-00	563-01110-00	564-01100-00	564-01110-00
G3 Vertical	563-01300-00	563-01310-00	564-01300-00	564-01310-00
G3 XL Vertical	563-01305-00	563-01315-00	564-01305-00	564-01315-00
G5 Dry Ice Condenser	563-01500-00	563-01510-00	564-01500-00	564-01510-00
G6 Vertical for Reflux	563-01600-00	563-01610-00	564-01600-00	564-01610-00

➤ Glassware

- Glassware set includes one 1,000-ml evaporating flask and one 1,000-ml receiving flask
- All glassware sets are also available with transparent plastic coating for added safety ¹
- Evaporating flasks and vapor tubes which come with a standard NS 29/32 joint size are also available with NS 24/29 ground joints
- Also available upon request: glassware sets for Soxhlet extraction and descending condenser system
- Innovative ventilation and replenishment valves without ground joints for easy handling and grease-free operation are included
- Glassware sets G3 and G6 allow for the automated distillation program AUTO_{accurate} on Hei-VAP Precision

XL Condenser


For maximum evaporation performance

The new XL condenser with industry-leading condensing surface of 2,200 cm² achieve up to 40 % higher evaporation rates. Available in G₃XL, G₃B XL and G₆² XL versions




Glassware set G₁

- Diagonal condenser for all standard distillations; the option with the **most affordable price**

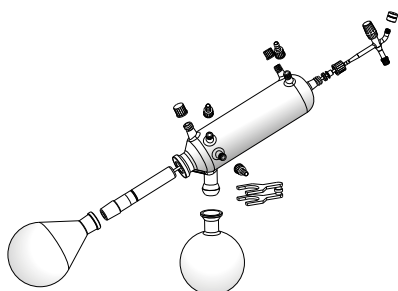


Glassware set G₃/G₃ XL

- Vertical condenser for all standard distillations; the option with the smallest space requirements, also available as XL version

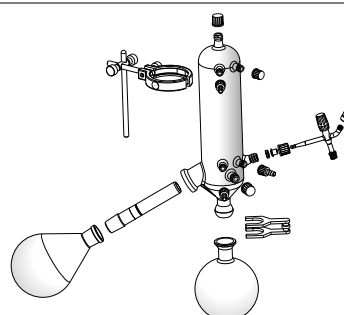


Glassware set G₅

- **Dry Ice** condenser for low-boiling solvents

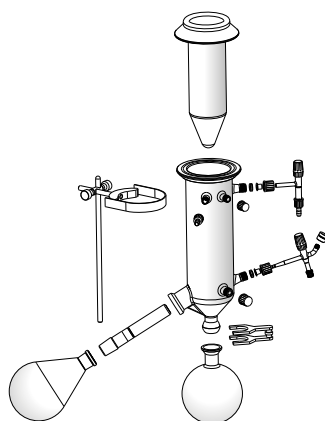

Glassware set G₆

- Vertical condenser; centerpiece with valve **for reflux distillation**

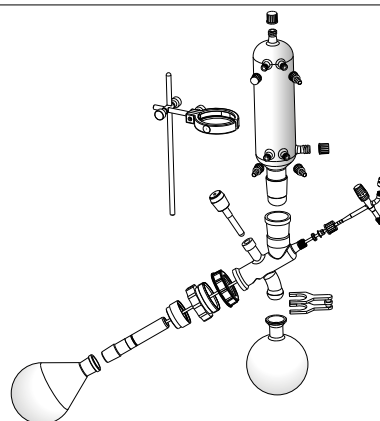

¹ Convenient for temperature range up to 80 °C

² for the usage of Deistimatic Bench-top

➤ Accessories



Evaporating flasks

	NS 29/32	NS 24/29
Evaporating flask* 50 ml	514-70000-00	514-70000-01
Evaporating flask 100 ml	514-71000-00	514-71000-01
Evaporating flask 250 ml	514-72000-00	514-72000-01
Evaporating flask 500 ml	514-73000-00	514-73000-01
Evaporating flask 1,000 ml	514-74000-00	514-74000-01
Evaporating flask 2,000 ml	514-75000-00	514-75000-01
Evaporating flask 3,000 ml	514-76300-00	514-76300-01
Evaporating flask 5,000 ml	514-78000-00	514-78000-03

* Evaporating flask NS 14.5 incl. adaptor (NS 29 or NS 24)




Receiving flasks

	standard	coated
Receiving flask - S 35/20 100 ml	514-81000-00	514-81000-02
Receiving flask - S 35/20 250 ml	514-82000-00	514-82000-02
Receiving flask - S 35/20 500 ml	514-83000-00	514-83000-02
Receiving flask - S 35/20 1,000 ml	514-84000-00	514-84000-02
Receiving flask - S 35/20 2,000 ml	514-85000-00	514-85000-02
Receiving flask - S 35/20 3,000 ml	514-87000-00	514-87000-02



Powder flasks

	NS 29/32	NS 24/29
Powder flask 500 ml	514-73200-00	514-73200-01
Powder flask 1,000 ml	514-74200-00	514-74200-01
Powder flask 2,000 ml	514-75300-00	514-75300-01

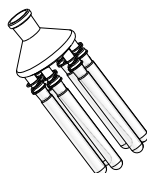

Spider

	NS 29/32	NS 24/29
Spider with 6 distillation sleeves NS 14,5 20 ml	15-300-005-09	15-300-005-12
Spider with 12 distillation sleeves NS 14,5 20 ml	15-300-005-10	15-300-005-13
Spider with 20 distillation sleeves NS 14,5 20 ml	15-300-005-11	15-300-005-06



Spider

	NS 29/32	NS 24/29
Spider with 5 flasks NS 24 50 ml	15-300-005-04	15-300-005-08
Spider with 5 flasks NS 24 100 ml	15-300-005-07	15-300-005-05



Evaporating cylinder

	NS 29/32	NS 24/29
Evaporating cylinder 500 ml	514-00070-00	514-00070-01
Evaporating cylinder 1,500 ml	514-00071-00	514-00071-01



Foam brake

	NS 29/32	NS 24/29
Foam brake with flask clamp	514-00004-00	514-00005-00

➤ Accessories



Guard hood

P/N 569-00010-00


- Strong and transparent guard hood to protect you from implosions
- Provides easy and immediate access to flask and bath



Guard shield

P/N 569-00020-00

- Strong and transparent guard shield for your protection while changing flasks
- Attached at the edge of the bath

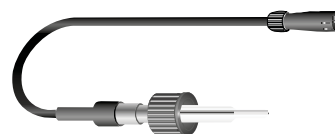

Vapor temperature sensor

P/N 569-00030-00

- For Hei-VAP Advantage and Hei-VAP Precision



AUTOaccurate-sensor

P/N 569-00040-00

- For Hei-VAP Precision
- Only in combination with glassware set G3 and G6



Tubing

P/N 591-35000-00 (package 12 m)

P/N 591-33000-01 (by meter)

- Tubing for vacuum and water connection. Available as package (12 meter piece incl. hose clamps) or in individual length (by meter)
- Inner diameter: 8 mm
- Outer diameter: 13 mm



Vacuum seal

P/N 23-30-01-01-30 (PTFE/FKM)

P/N 23-30-01-06-70 (PTFE/FFKM)

- Standard PTFE / FKM spare seal
- Ultra resistant PTFE / FFKM seal for applications with highly abrasive solutions



Aqua-Stop

P/N 569-00500-00

The integrated valve closes automatically by turning off the evaporator. The system saves money by reducing your consumption of cooling tap water. Moreover the safety mechanism avoids flooding in your laboratory. (Only for Hei-VAP Advantage and Hei-VAP Precision)



Spare clamping sleeve

P/N 23-30-01-05-31

Spare vapor tube NS 29/32

P/N 514-00000-01



Heating bath liquid

P/N 569-00600-00

- Heating bath liquid up to 220 °C (5 kg)

➤ Vacuum Pumps

Advantages of all chemical-resistant diaphragm vacuum pumps



- All components which come in contact with media are built from chemically resistant fluoropolymer
- The durable PTFE design guarantees an superior diaphragm lifespan
- The head cover and clamping disc have a stable core made of metal which offers unsurpassed long-term performance of your operational parameters
- The direct pump drive (without belt) is exceptionally quiet, creates a very low vibrational environment and reduces the need of wear and tear parts to a minimum
- The gas ballast valve has been optimized to prevent media condensing in the pump


➤ RPM-Regulated Vacuum Pumps

➤ **Highest distillation rates possible and reducing process time by 30 %**

➤ **Stable distillation conditions account for solvent recovery up to 99 %**

YOUR ADVANTAGES

- Speed adjustment stops completely if vacuum is reached and runs intermittently at very low speed to hold vacuum constant
- Increased diaphragm performance life due to periods of non-operation
- Low energy consumption due to intermittent operation
- Significant noise and vibration reduction due to periods of non-operation
- Operates without vacuum valve
- Most precise vacuum level at all times

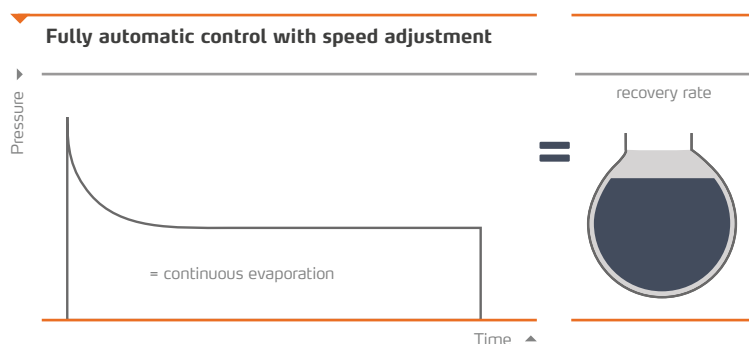

RPM regulated pumps adjust their speed individually when the required vacuum level is reached

➤ Valve-Regulated Vacuum Pumps

Standard vacuum pumps can be controlled manually or via valve operated vacuum controllers

Rotavac Valve Tec



Rotavac Valve Tec

P/N 591-00160-00

Condenser

P/N 591-00083-00

**For all
Hei-VAP models**

Diaphragm pump		Two-stage
Suction capacity	m ³ /h	0.75
Ultimate vacuum	mbar	12
Combinable with a condensor		yes
Power input	Watt	80
Weight	kg	6.0
Dimensions L/W/H (without condenser)	mm	145 / 315 / 169

Specifications per unit

Recommended for solvents with low or medium boiling points

Vacuum can be controlled via valve operated or manual vacuum controllers

Rotavac Valve Control



Rotavac Valve Control

P/N 591-00130-00

Condenser

P/N 591-00083-00

**For all
Hei-VAP models**

Diaphragm pump		Two-stage
Suction capacity	m ³ /h	2.0
Ultimate vacuum	mbar	7
Combinable with a condensor		yes
Power input	Watt	180
Weight	kg	12.8
Dimensions L/W/H (without condenser)	mm	195 / 245 / 310

Suction capacity for up to three rotary evaporators at the same time

Vacuum can be controlled via valve operated or manual vacuum controllers

Depending on your application you can switch on and off the vacuum pump via switchbox

➤ RPM-Regulated Vacuum Pumps

Controllable vacuum pumps for digital regulation of the vacuum via operation panel or directly at the unit (depending on version - see below)

Rotavac Vario Tec



Rotavac Vario Tec

P/N 591-00171-00

Condenser

P/N 591-00084-00

**For
Hei-VAP Precision**

Two-stage

1.0

12

yes

160

4.3

156 / 236 / 196

Recommended for solvents with low or medium boiling points

Rotavac Vario Control



Rotavac Vario Control

P/N 591-00141-00

Condenser

P/N 591-00084-00

**For
Hei-VAP Precision**

Three-stage

1.7

2

yes

160

5.4

167 / 236 / 196

Even if the gas ballast valve is open an excellent ultimate vacuum is reached when working with easily condensable vapors. This makes distilling high boiling point solvents such as DMF or DMSO possible at low bath temperatures

Rotavac Vario Pumping Unit



Rotavac Vario Pumping Unit

P/N 591-00142-00

Condenser

P/N 591-00084-00

Fully controllable stand-alone pumping unit including vacuum controller. For Hei-VAP Value and Hei-VAP Advantage

Three-stage

1.7

5

yes

160

6.0

193 / 263 / 299

Precise vacuum control avoids bumping and eliminates foaming of your evaporation solution

Due to automatic vacuum supply to the process parameters you have more time for important laboratory activities

➤ Vacuum Accessories

Switchbox

Connect up to three rotary evaporators to one vacuum pump



Switchbox

P/N 569-00400-00

- Connect up to to three rotary evaporators to one pump and control vacuum individually
- Pump switches off automatically at required vacuum
- Reduces energy consumption of pump
- Increases performance life of the pump due to a complete switch-off if vacuum is reached
- Vacuum valve and non-return valve are required for every evaporator if you work with different pressures at the same time
- Delivery includes three non-return valves
- Designed for Hei-VAP Precision
- Weight: 0.6 kg
- Dimensions: L80 / W160 / H45 mm

Manual vacuum controller

For Hei-VAP Value and Hei-VAP Advantage



Manual vacuum controller

P/N 591-26000-00

- Control your house vacuum or any other vacuum source
- Pressure range from 0 to 1,020 mbar
- Scaling on display in steps of 50 mbar
- Connection for tubes with 8 mm inner diameter
- Support rod included
- A vacuum valve is not required
- Dimensions: L80 / W80 / H150 mm

Vacuum valve

Required to control vacuum for Hei-VAP Precision models with valve-regulated pumps



Vacuum valve

P/N 569-00060-00

- External attachment design for easy cleaning
- Can be combined with the Woulff Bottle to one unit
- Not required if you utilize a Manual vacuum controller

Woulff bottle

P/N 569-00070-00

- Additional recovery of solvent to protect vacuum valve and vacuum pump
- Compatible with all Hei-VAP configurations
- Direct connection to the vacuum valve is possible
- Features connection to vacuum controller
- Volume of 250 ml

Chillers

RotaChill Large

Highest cooling capacity



RotaChill Large

P/N 591-00240-00

- Highest cooling capacity for any application
- Temperature range from -20 °C to +40 °C
- Temperature control accuracy of ± 0.1 °C
- Digital display for temperature setting and reading out actual temperature
- Cooling capacity at +20 °C: 1,290 W
- Footprint of chiller: L607 / W254 / H483 mm

Rotacool Mini

Ideal for standard applications



Rotacool Mini

P/N 591-00210-00

- Chiller for all standard applications
- Temperature range from -10 °C to +40 °C
- Temperature control accuracy of ± 1 °C
- Digital display for temperature settings and reading out actual temperature
- Cooling capacity at +20 °C: 580 W
- Footprint of chiller: L360 / W310 / H520 mm

RotaChill Small

The affordable model



RotaChill Small

P/N 591-00220-00

- Temperature range from 0 °C to +100 °C
- Precise temperature control accuracy of ± 0.2 °C
- The RotaChill Small should be placed directly besides the evaporator to ensure a sufficient pump performance
- Footprint of chiller: L450 / W210 / H570 mm

		RotaChill Large	Rotacool Mini	RotaChill Small
Temperature range	[°C]	-20 to +40	-10 to +40	0 to +100
Coolant volume	[L]	2.7	3	6
Cooling capacity* at: +20 °C	[W]	1,290	580	200
0 °C		750	370	140
-10 °C		-	80	-
-20 °C		260	-	-

* at 20 °C ambient temperature

➤ Equipment Combinations

With the Hei-VAP series you can easily assemble your individual laboratory equipment - exactly to your requirements.

Please refer to the overview by Hei-VAP units for the possible configurations or contact us - we are happy to be at your service and will advise you in detail.

➤ Example configuration

Rotary Evaporator



Hei-VAP Advantage

Vacuum source



Rotavac Vario
Pumping Unit



Accessories



Vacuum controller

Vapor temperature sensor
(all glassware sets)

Vacuum valve

included



included in pumping unit, no further configuration required

optional



optional available as accessory

not required



not required due to rpm regulation

As an aid to decision-making we have created different packages for you, from basic equipment up to automatic distillation with RPM-regulated vacuum pumps - for details see the following pages.

➤ Overview of Hei-VAP units

Rotary Evaporator



Hei-VAP Value / Hei-VAP Value Digital

Vacuum source

Rotavac Vario
Pumping Unit

included ➤

not required ➤

Rotavac Valve Control
Rotavac Valve Tec
House vacuum

optional ➤

Accessories

Vacuum controller

Vacuum valve

Manual vacuum
controller



Hei-VAP Advantage

Rotavac Vario
Pumping Unit

included ➤

optional ➤

not required ➤

Rotavac Valve Control
Rotavac Valve Tec
House vacuum

optional ➤

Vacuum controller

Vapor temperature sensor
(all glassware sets)

Vacuum valve

Manual vacuum
controller

Vapor temperature sensor
(all glassware sets)



Hei-VAP Precision

Rotavac Valve Control
Rotavac Valve Tec
House vacuum

included ➤

required ➤

optional ➤

Rotavac Vario Control
Rotavac Vario Tec

included ➤

not required ➤

optional ➤

Vacuum controller

Vacuum valve

Auto_{accurate}-sensor
(glassware G3/G6 only)

Vapor temperature sensor
(all glassware sets)

Vacuum controller

Vacuum valve

Auto_{accurate}-sensor
(glassware G3/G6 only)

Vapor temperature sensor
(all glassware sets)



Rotavac Valve Control
with
condenser



Rotavac Valve Tec
with
condenser



Rotavac Vario Control
with
condenser



Rotavac Vario Tec
with
condenser



Rotavac Vario
Pumping Unit
with condenser



Manual vacuum
controller



Vacuum valve
with Woulff bottle



Auto_{accurate}-
sensor



Vapor
temperature sensor

➤ Hei-VAP Packages



Hei-VAP Silver 4 Package

P/N 560-91302-00



Hei-VAP Gold 3 Package

P/N 562-91307-00

Hei-VAP Silver Packages

Standard equipment for your lab

Hei-VAP Silver 4 Package

P/N 560-91302-00

Hei-VAP Value Digital with hand lift and vertical glassware set G3

- Vacuum pump Rotovac Valve Tec
- Manual vacuum controller
- Tube set

Hei-VAP Silver 5 Package

P/N 560-81312-00

Hei-VAP Value Digital with hand lift and coated vertical glassware set G3B

- Vacuum pump Rotovac Valve Control
- Manual vacuum controller
- Tube set

Hei-VAP Silver 6 Package

P/N 560-91309-20

Hei-VAP Value Digital with hand lift and vertical glassware set G3

- Vacuum pump Rotovac Valve Tec
- Manual vacuum controller
- Chiller RotaChill Small
- Tube set

Hei-VAP Gold Packages

For standard distillations with vacuum control

Hei-VAP Gold 1 Package

P/N 561-41300-00

- Hei-VAP Advantage with hand lift and vertical glassware set G3
- Vapor temperature sensor
- Vacuum pump Rotovac Valve Tec
- Manual vacuum controller
- Tube set

Hei-VAP Gold 2 Package

P/N 562-11310-00

- Hei-VAP Advantage with motor lift and coated vertical glassware set G3B
- Vapor temperature sensor
- Vacuum pump Rotovac Vario Pumping Unit
- Tube set

Hei-VAP Gold 3 Package

P/N 562-91307-00

- Hei-VAP Advantage with motor lift and vertical glassware set G3
- Vacuum pump Rotovac Valve Tec
- Manual vacuum controller
- Chiller Rotacool Mini
- Tube set

Hei-VAP Gold 4 Package

P/N 564-61310-00

- Hei-VAP Precision with motor lift and vertical glassware set G3B
- Vacuum pump Rotovac Valve Tec with condenser
- Vacuum valve
- Tube set



Hei-VAP Platinum 2 Package
P/N 563-51310-00



Hei-VAP Platinum 5 Package
P/N 564-31318-00

Hei-VAP Platinum Packages

Automatic distillations with Hei-VAP Precision hand lift models and valve-regulated vacuum pumps

Hei-VAP Platinum 1 Package

P/N 563-61300-00

- Hei-VAP Precision with hand lift and vertical glassware set G3
- Vapor temperature sensor
- Vacuum valve
- Vacuum pump Rotovac Valve Tec with condenser
- Tube set

Hei-VAP Platinum 2 Package

P/N 563-51310-00

- Hei-VAP Precision with hand lift and coated vertical glassware set G3B
- AUTOaccurate-sensor
- Vacuum valve
- Vacuum pump Rotovac Valve Control with condenser
- Tube set

Hei-VAP Safety Upgrade

For highest safety

Hei-VAP Safety Upgrade

P/N 569-00311-00

Additional guard hood and coating of glassware for your Hei-VAP Package



Hei-VAP Platinum Packages

Automatic distillations with Hei-VAP Precision motor lift models and RPM-regulated vacuum pumps

Hei-VAP Platinum 3 Package

P/N 564-21300-00

- Hei-VAP Precision with motor lift and vertical glassware set G3
- AUTOaccurate-sensor
- Vacuum pump Rotovac Vario Tec with condenser
- Woulff bottle
- Tube set

Hei-VAP Platinum 4 Package

P/N 564-31310-00

- Hei-VAP Precision with motor lift and coated vertical glassware set G3B
- AUTOaccurate-sensor
- Vacuum pump Rotovac Vario Control with condenser
- Woulff bottle
- Tube set

Hei-VAP Platinum 5 Package

P/N 564-31318-00

- Hei-VAP Precision with motor lift and coated vertical glassware set G3B
- AUTOaccurate-sensor
- Vacuum pump Rotovac Vario Control with condenser
- Woulff bottle
- Chiller RotaChill Large
- Tube set

Hei-VAP Platinum 5 XL Package

P/N 564-31319-00

- Hei-VAP Precision with motor lift and coated vertical glassware set G3B XL
- AUTOaccurate-sensor
- Vacuum pump Rotovac Vario Control with condenser
- Woulff bottle
- Chiller RotaChill Large
- Tube set

Technical Specifications - Rotary Evaporators

Model	Hei-VAP Value (HL) / Value Digital (HL)	Hei-VAP Advantage (HL)	Hei-VAP Advantage (ML)	Hei-VAP Precision (HL)	Hei-VAP Precision (ML)
Hand / Motor lift	hand lift	hand lift	motor lift	hand lift	motor lift
Height adjustment speed	manual	manual	30 mm/s	manual	30 mm/s
Height adjustment (mm)	155	155	155	155	155
Rotation speed (rpm)	10 – 280	10 – 280	10 – 280	10 – 280	10 – 280
Rotation speed setting	scale	3.5" LCD display	3.5" LCD display	4.3" LCD color graphic display	4.3" LCD color graphic display
Drive	brushless DC motor with electronic speed control				
Heating capacity (W)	1,300	1,300	1,300	1,300	1,300
Temp. range heating bath (°C)	20 – 210	20 – 210	20 – 210	20 – 210	20 – 210
Temperature accuracy bath (°C)	±1	±1	±1	±1	±1
Overheat protection bath (°C)	cut-off at 5 °C over set temperature via separate Pt 1000				
Bath temperature setting***	scale / digital	3.5" LCD display	3.5" LCD display	4.3" LCD color graphic display	4.3" LCD color graphic display
Heating bath temperature control	electronic / digital	electronic / digital	electronic / digital	electronic / digital	electronic / digital
Secondary overtemp. cut-off (°C)	250	250	250	250	250
Material heating bath	stainless steel AISI 316L	stainless steel AISI 316L	stainless steel AISI 316L	stainless steel AISI 316L	stainless steel AISI 316L
Diameter heating bath (mm)	255	255	255	255	255
Volume heating bath (L)	4.5	4.5	4.5	4.5	4.5
Display vapor temperature required accessories	no	yes vapor temperature sensor	yes vapor temperature sensor	yes vapor temperature sensor	yes vapor temperature sensor
Integrated vacuum controller	no	no	no	yes	yes
Timer	no	yes	yes	yes	yes
Vacuum program SET _{pressure} **	no	no	no	yes	yes
Vacuum program AUTO _{easy} **	no	no	no	yes	yes
Vacuum program AUTO _{accurate} ** requires AUTO _{accurate} -sensor	no	no	no	yes	yes
Remote control function					
- Lift	no	no	yes	no	yes
- Start rotation	yes	yes	yes	yes	yes
- Start vacuum	no	no	no	yes	yes
- Start heating bath	yes	yes	yes	yes	yes
Programmable ramps	no	no	no	yes	yes
Supply power (W)	1,400	1,400	1,400	1,400	1,400
Weight (w/o glassware assembly) (kg)	16	16	16	17	17
Dimensions with glassware set G3 (mm) Width on extended heating bath (l x w x h)	477 x max. 739 x 845	477 x max. 739 x 845	407 x max. 739 x 845	523 x max. 739 x 845	453 x max. 739 x 845
Dimensions with glassware set G3 XL (mm) Width on extended heating bath (l x w x h)	477 x max. 739 x 920	477 x max. 739 x 920	407 x max. 739 x 920	523 x max. 739 x 920	453 x max. 739 x 920
Protection class*** (DIN EN 60529)	IP 20	IP 20	IP 20	IP 20	IP 20
Permissible ambient conditions	5 – 31 °C at 80 % rel. humidity 32 – 40 °C decreasing linearly up to max. 50 % rel. humidity	5 – 31 °C at 80 % rel. humidity 32 – 40 °C decreasing linearly up to max. 50 % rel. humidity	5 – 31 °C at 80 % rel. humidity 32 – 40 °C decreasing linearly up to max. 50 % rel. humidity	5 – 31 °C at 80 % rel. humidity 32 – 40 °C decreasing linearly up to max. 50 % rel. humidity	5 – 31 °C at 80 % rel. humidity 32 – 40 °C decreasing linearly up to max. 50 % rel. humidity
Condensing surface (cm ²)	1,400	1,400	1,400	1,400	1,400
Condensing surface XL Condenser (cm ²)	2,200	2,200	2,200	2,200	2,200

Rate of evaporation* (L/h)	Toluene	Acetone	Ethanol	Isopropanol (IPA)	Water
G3 Vertical	6.5	6.0	2.5	4.0	1.4
G3 XL Vertical	8.5	7.5	3.5	4.9	1.4
Improvement with XL condenser (%)	31	25	40	23	system limit reached with 1.4 L/h

* With maximum condenser efficiency
 ** Only in combination with vacuum systems

*** IP 67 cable connection on heating bath
 **** Hei-VAP Value Digital: digital display

Standard supply voltage: 230 V - other voltages upon request, please specify for order

Certificate

**To confirm the ability for
continuous operation
of the Hei-VAP series Rotary Evaporators**

Hei-VAP series Rotary Evaporators feature overtemperature safety circuits according to DIN EN 61010-1:2001 and DIN EN 61010-2-010:2014 and therefore is designed for continuous operation.

This statement is made under the precondition that all units are operated in accordance with the operation manual and in accordance with good practice standards for safety in laboratories, rules for accident preventions, and compliance with directions on hazardous materials.

S c h w a b a c h , J a n u a r y 2 0 1 8



Stefan Peters
Research and Development Manager



Marcell Sarré
Quality Manager